

Please Note:

- ⚠ Do not look directly at the laser light as this may cause harm to your eyes. It is a normal phenomenon that your scalp becomes warm during usage of the device.
- ⚠ WARNING: Do not modify this equipment without authorization of the manufacturer.

04 Main Technical Specification

Product name	Laser Cap 272
Laser diodes	272 pcs * 5 mW(each)
Laser wavelength	650nm
Working mode	4 zoned irradiation mode: F (Forehead), C (Crown), L (Lower Back), or H (Full Scalp)
Output power	5 mW each
Treatment time	15 minute (continuous mode) / 30 minute (pulse mode)
Battery capacity	5200mAh



Laser Hair Growth System

- ◆ 272 medical-grade lasers
- ◆ 1360 mW output power
- ◆ Targeted zoned irradiation



Contact us

Manufacturer: Hubei YJT Intelligent Technology Group Co., Ltd.

ADD: 8F,Block7,No.62 Optical Valley Road,Wuhan,China.

Tel: +86-27-51231823

Email: weldon@yjt.com.cn

Site: www.yjtlascure.com
www.yjthnc.com



01 Description


Laser cap 272 is indicated to promote hair growth in males with androgenic alopecia who have Norwood-Hamilton classifications of IIa-V or females with androgenic alopecia who have Ludwig-Savin Classifications of I-III and both with Skin Phototypes I-IV.

Laser 650nm can penetrate the surface of the activities of organisms and then increase their ATP. It can also promote the tissue of regeneration, speed up the dissipation of fatigue- and promote the expansion of capillaries. With its light therapy, it can improve the absorbency of hair follicles, so they can easily absorb nutrients, and therefore activate cells, repair damaged roots, and prompt regeneration of the hairs.

Laser cap 272 is a dome-shaped low level laser therapy (LLLT) device is designed to promote hair growth in women and men by exposing the entire scalp to the photobiostimulation of 272 pcs of 650nm laser diodes. The Helmet is designed with an outer plastic cover and a protective inner liner (containing the electronics and laser array).

02 Indications

- Promotes Hair Growth
- Reduces Shedding
- Strengthens Hair
- Reverses the Thinning Process
- Increases Density & Fullness
- Energizes & Awakens Follicles
- Rejuvenates Aging Hair
- Restores the Hair Growth Cycle
- Revitalizes Damaged Hair
- Frontal patterns of hair loss
- Norwood-Hamilton Classification IIa to V

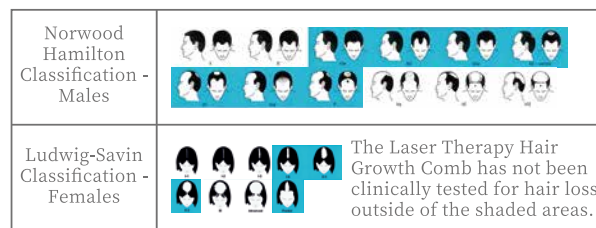

Before

After



Treating time :

The Laser Hair Growth System should be used 3 times per week, 15 minutes or 30 minutes per session, on non-consecutive days (e.g., Monday, Wednesday, Friday). Use for 17 weeks for initial results. After that, use as needed for maintenance, not exceeding 30 minutes per day. Select one mode per session.

Package Content

When you open the gift box of Laser Hair Growth System, you'll find it contains the following components:

- ◆ Inner laser Cap *1
- ◆ Controller *1
- ◆ User Manual *1
- ◆ Outer Cloth Cap *1
- ◆ Power Adapter *1
- ◆ Colorful Box *1

03 Operation

- 1) Remove all components from the gift box;
 - 2) Make sure the Cap is clean, dry and intact before beginning the treatment. You can clean the Cap with a soft, slightly moistened cloth. This is done in order to remove hairs and dust from the Cap;
 - 3) Short press the on/off button until the controller lights up, then press the "T" button to select your session: 15-minute continuous mode (CW) or 30-minute pulse mode (PW).
- The cap features targeted zonal irradiation. Press the "M" button to adjust the treatment area (4 modes available): F (Forehead), C (Crown), L (Lower Back), or H (Full Scalp)
- 4) During usage, Power display Light three lights, current power 80%-100%; Light two lights, when the battery is 50%-80%; Light a light, the current power is less than 50%.
 - 5) You may clean the surface of the Cap with a medicinal alcohol wipe or slightly moistened cloth after each session.
 - 6) The Cap has a built-in battery. When the power is less than 20%, the indicator flashes, and has a drip prompt sound, and need to charge.

INTELLIGENT ZONAL IRRADIATION

4 Zonal Irradiation Modes: Forehead, Crown, Back of the Head, Full Scalp

